


# Technical Specifications

## MEGA-PURE MP-11A

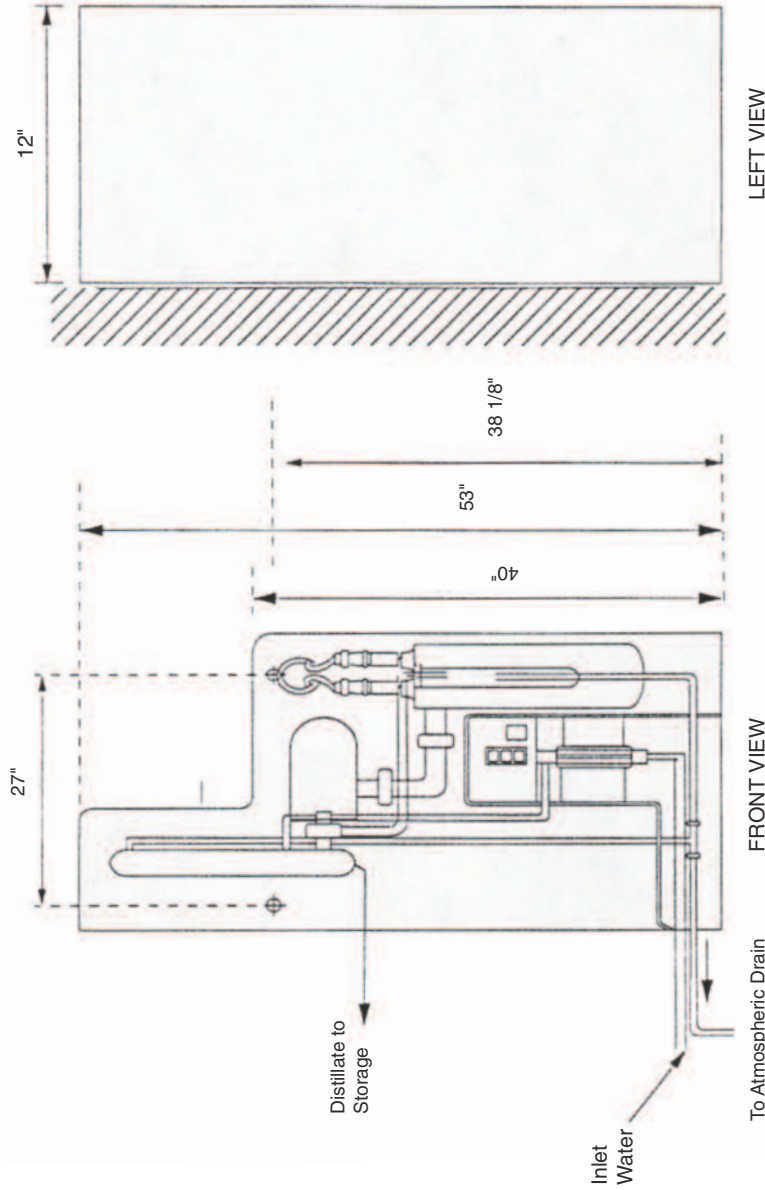

**SPECIFICATIONS**  
 Distilled Water Flow - 13 liters/hr.  
 Feedwater Flow - 130 liters/hr.  
 Feedwater Pressure - 20-100 psig  
 Operating Weight - 125 lbs.

**ELECTRICAL**  
 A440118 - 240V, 1 phase, 42 amp, 50/60 Hz, 10 KW  
 A440117 - 208V, 1 phase, 49 amp, 50/60 Hz, 10 KW

- NOTES:**
- CAUTION:** Refer to appropriate Barnstead Equipment manuals before installing, operating or servicing this equipment.
  - CAUTION:** Water from this unit is for a specific use. For best results, carefully follow the instructions in the appropriate Barnstead manuals.
  - All service piping and wiring by customer to be in accordance with national and local codes.
  - Installer to attach unit to wall using suitable fasteners not supplied by Barnstead.
  - Do not mount electrical equipment or outlets below this unit.
  - For additional information see owners manual.
  - A 6" clearance is required around the perimeter of the still.
  - Unit requires an atmospheric drain.
  - Unit is supplied with a cord and plug for a Hubbel #9460 receptacle.

<b>Drawing No.</b>	RS 6003-93 (1741.4-92)
<b>Description</b>	MEGA-PURE - MP11A
<b>Scale</b>	Not to Scale
<b>Date</b>	05/18/04

**Barnstead International • 2555 Kerper Blvd., P.O. Box 797 • Dubuque, Iowa 52004-0797 • 1-800-553-0039**