


Technical Specifications

MEGA-PURE MP-6A



LEFT VIEW

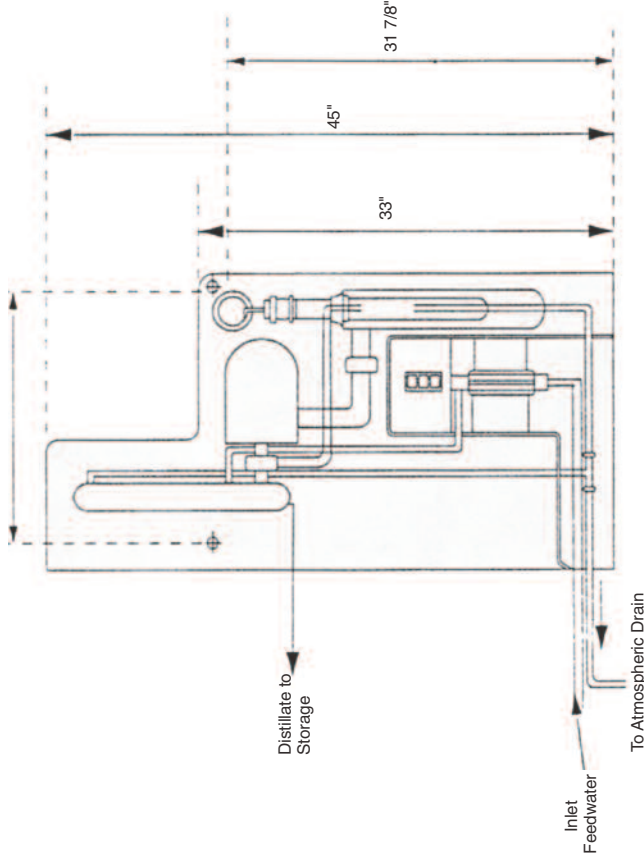

FRONT VIEW

SPECIFICATIONS

Distilled Water Flow - 6 liters/hr.
 Feedwater Flow (tap feed) - 60 liters/hr.
 Feedwater Pressure - 20-100 psig
 Operating Weight - 110 lbs.

ELECTRICAL

A440518 - 240V, 1 phase, 21 amp, 50/60 Hz, 5 KW
 A440697 - 208V, 1 phase, 25 amp, 50/60 Hz, 5 KW

NOTES:

- CAUTION:** Refer to appropriate Barnstead Equipment manuals before installing, operating or servicing this equipment.
- CAUTION:** Water from this unit is for a specific use. For best results, carefully follow the instructions in the appropriate Barnstead manuals.
- All service piping and wiring by customer to be in accordance with national and local codes.
- Installer to attach unit to wall using suitable fasteners (customer supplied).
- Do not mount electrical equipment or outlets below this unit.
- For additional information see owners manual.
- A 6" clearance is required around the perimeter of the still.
- Unit requires an atmospheric drain.
- Unit is supplied with a cord and plug for a Hubbel #5661 receptacle.

Drawing No.	RS 6002-93 (1741.2-92)
Description	MEGA-PURE - MP6A
Scale	Not to Scale
Date	02/26/04

Barnstead International • 2555 Kerper Blvd., P.O. Box 797 • Dubuque, Iowa 52004-0797 • 1-800-553-0039