


3/4" male NPT fitting with pressure and potable tap

- 1 Double Pretreatment for Barnstead Pacific TII/RO water systems (Item No. 09.4.000)
- 2 Barnstead Pacific TII/RO pure water system
- 3 Barnstead 100L Reservoir for Pacific TII/RO water systems
- 4 OPTIONAL Hand Dispenser can be utilized up to 3m (108") from water system
- 5 Holes are for optional wall mounting of water system with supplied hardware


Legend:

Notes:

1 Dual dimensions are in metric and inch

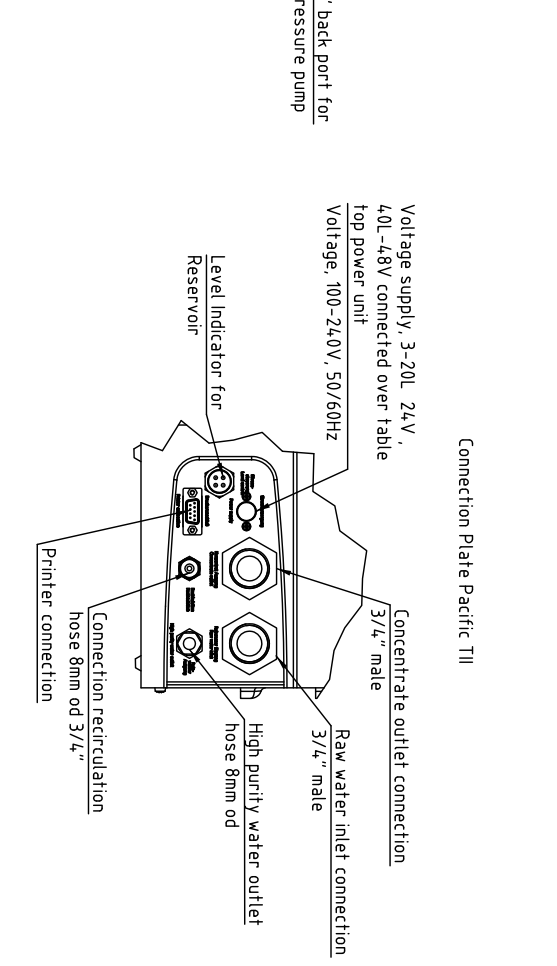

Connections from Tank

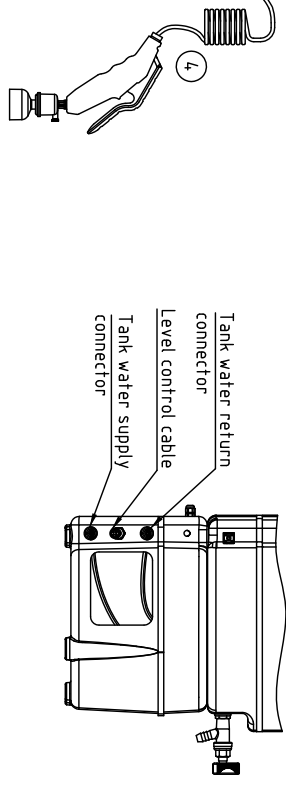

ELECTRICAL REQUIREMENTS		100-240 V, 50/60 Hz, 2-1A, up to 5 ft from unit
WATER CONNECTIONS		3/4" NPT with manual shut off valve recommended
MIN/MAX INLET PRESSURE		15 - 85 PSI (0.1 - 6 bar)
RECOMMENDED FEED TEMPERATURE		2 - 35°C
RECOMMENDED FEED WATER TYPE		Tap/Potable Water
DRAIN		An atmospheric drain must be available within 5 feet of the final mounting location

Client / Town:		Project:		Scale	
	Date	Name			
	Drawn	11.09.2014	Avadjeva		
	Checked	11.09.2014	ODP		
Dimension sheet		Pretreatment, Pacific, 100L Tank			
MB119		Sheet			
		Sheets			

c	Vorbehandlung	26.06.2015	M.A
b	geändert	26.11.2014	M.A
a	geändert	14.11.2014	M.A
	Amendment	Date	Name

Thermo
SCIENTIFIC

Thermo Electron LED GmbH
Stübeln 3, Düsseldorf
Tel. +49 (0)21 06 99-0